Target 3D

Instructions for import of the implant library
Target 3D in 3Shape Software



3shape 

How to import the Target 3D Library


Library is imported through the *Control Panel* of your 3Shape system

1. Choose *Import/Export*

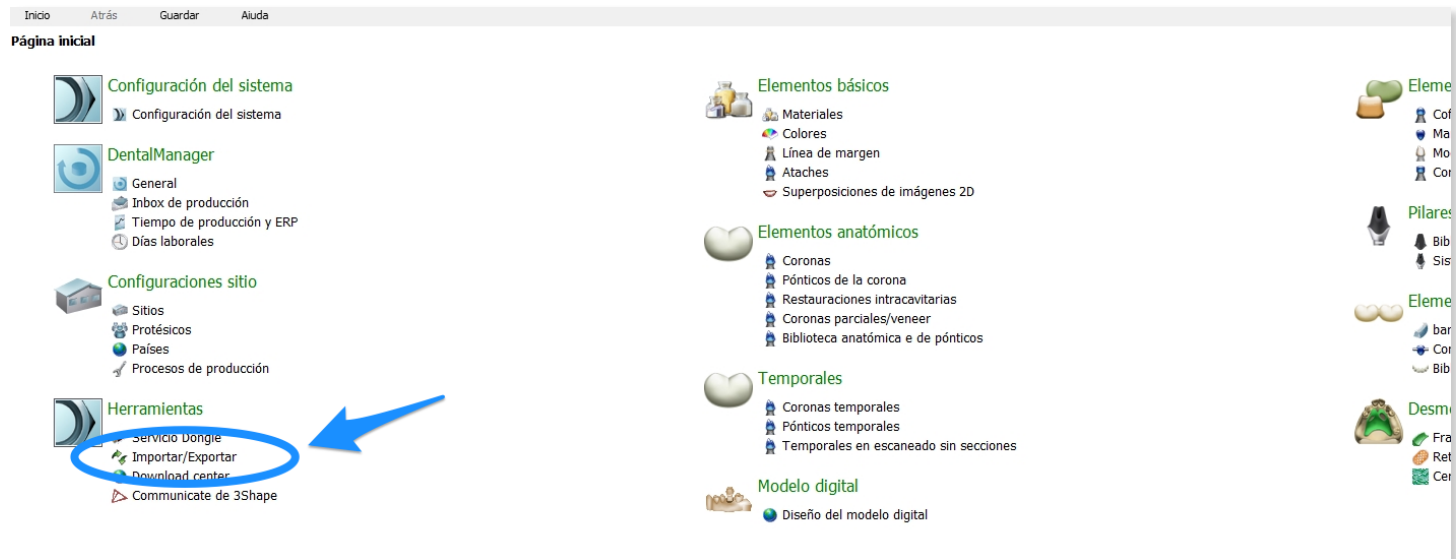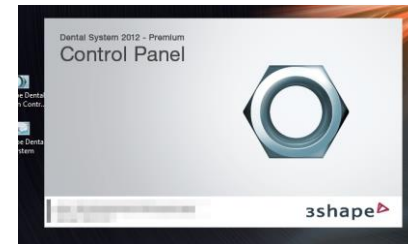
How to import the Target 3D Library

2. Choose *import material*

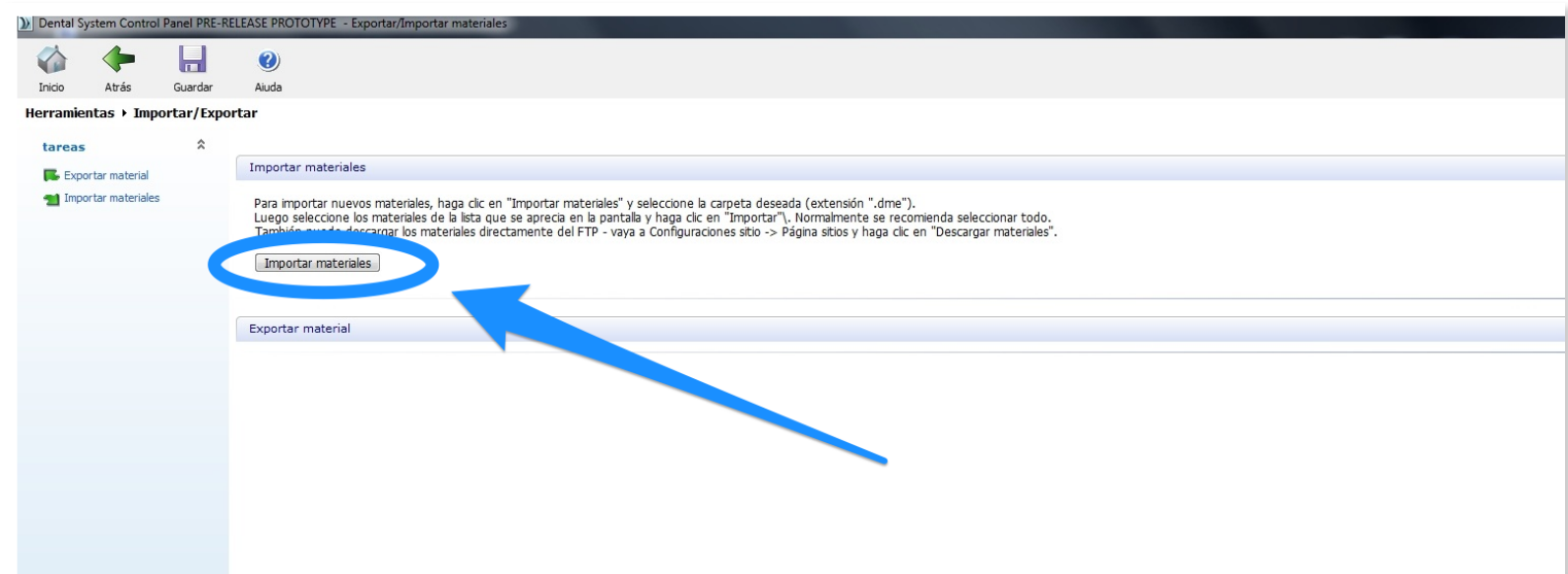
How to import the Target 3D Library

3. Choose and open Target 3D library

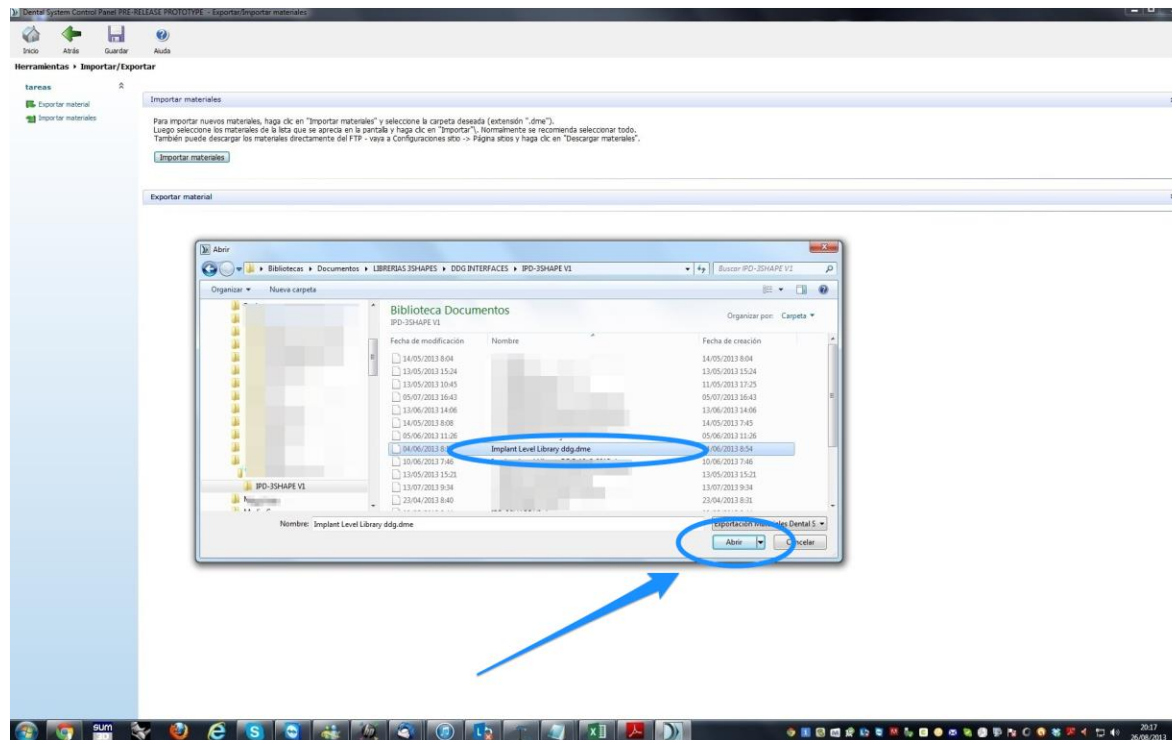
How to import the Target 3D Library

4. Import the Target 3D library

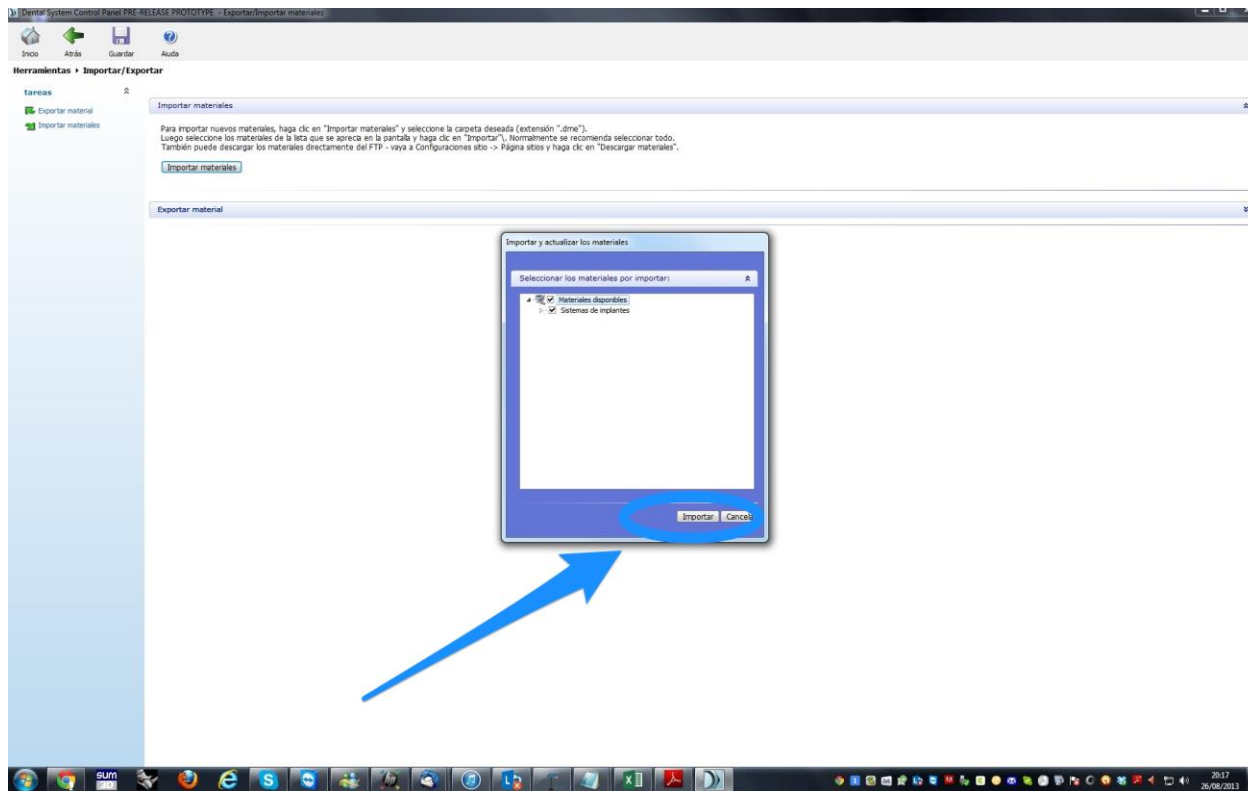How to import the Target 3D Library

5. Confirm succesful import

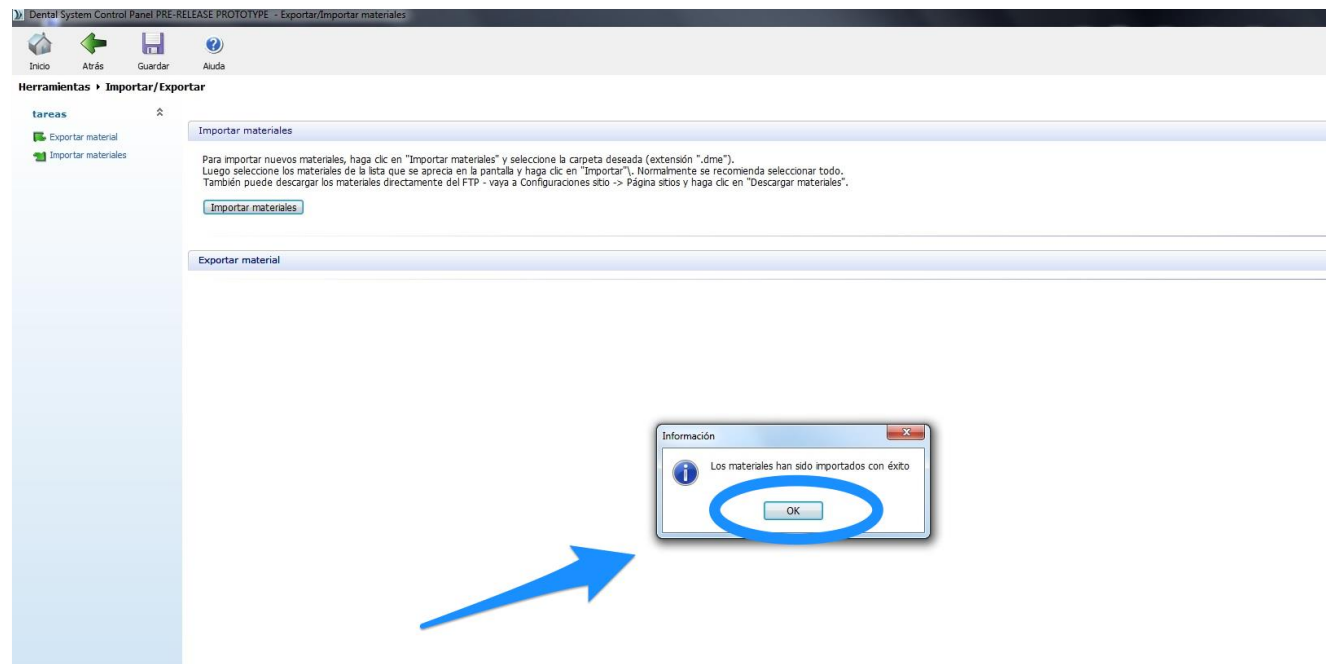How to import the Target 3D Library

6. Save Target 3D and confirm

The screenshot displays the 'Dental System Control Panel PRE-RELEASE PROTOTYPE - Exportar/Importar materiales'. The interface includes a top navigation bar with 'Inicio', 'Atrás', 'Guardar', and 'Ayuda' buttons. A blue circle highlights the 'Guardar' button. Below this, the 'Herramientas > Importar/Exportar' section is visible. The 'Importar materiales' button is highlighted with a blue arrow. The main content area contains instructions for importing materials and a list of tasks. A second screenshot shows the 'Guardar' button in the top navigation bar highlighted with a blue circle and arrow. Below it, a 'Confirmar guardar' dialog box is shown with the 'SI' button highlighted by a blue circle and arrow.

Dental System Control Panel PRE-RELEASE PROTOTYPE - Exportar/Importar materiales

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Importar materiales

Exportar material

Dental System Control Panel PRE-RELEASE PROTOTYPE

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Confirmar guardar

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
SI No Cancelar

Display of the imported Target 3D Library in the 3Shape Manager

That's all. Easy.

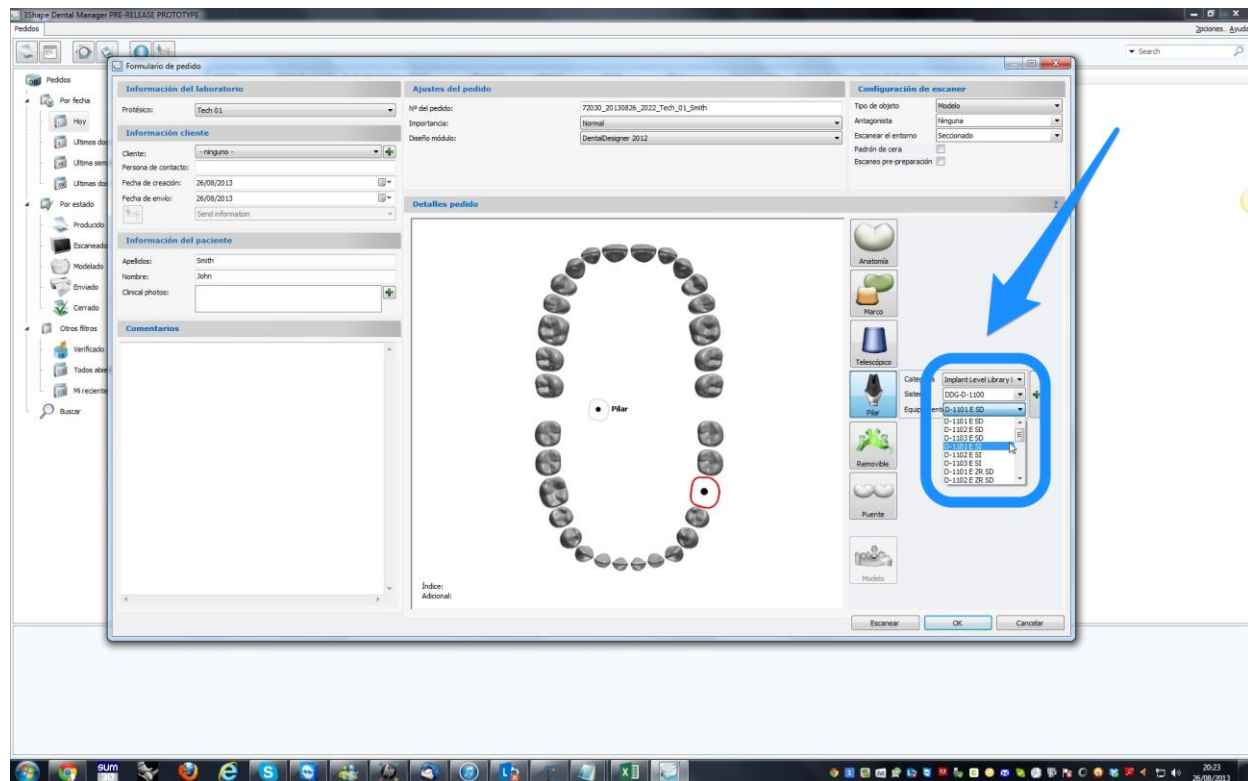Target 3D

3shape

Thanks for Watching

Target 3D Support

Target 3D